

SCALING



Used to remove light and medium supragingival calculus and bacterial plaque.

G1 GS1 GN1
GD1 GM1



Flat edge is used to remove heavy supragingival calculus.

G2 GS2 GN2
GD2



Used to remove calculus and bacterial plaque from supragingival interdental areas.

G3 GS3 GN3
GD3 GM3




Used to remove all supragingival calculus and bacterial plaque.

G4 GS4 GN4
GD4 GM4



Used to remove all supragingival calculus and bacterial plaque.

G5 GS5 GN5
GD5 GM5

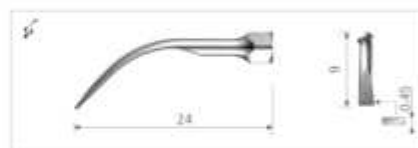

Used to remove heavy supragingival calculus.

G6 GS6 GN6
GD6 GM6

PERIODONTICS



Used to remove subgingival calculus.

P1 PS1 PN1
PD1 PM1




Curved to the left, for treating furcations, interdental calculus and very narrow spaces.

P2L PS2L PN2L
PD2L



Curved to the left, diamond coated, for treating furcations, interdental calculus and very narrow spaces.

P2LD PS2LD PN2LD
PD2LD

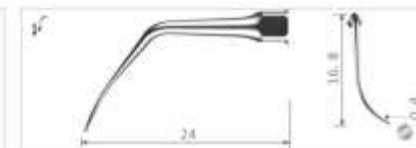

Curved to the right, for treating furcations, interdental calculus and very narrow spaces.

P2R PS2R PN2R
PD2R



Curved to the right, diamond coated, for treating furcations, interdental calculus and very narrow spaces.

P2RD PS2RD PN2RD
PD2RD



Used to remove subgingival calculus.

P3 PS3 PN3
PD3 PM3

EMS SAT SIR NSK



Diamond coated, used to smooth the surface of root during periodontal flap surgery.

P3D
PD3D



Used to remove subgingival calculus.

P4 PS4 PN4
PD4



Diamond coated; a powerful tool for the removal of calcification, this tip can be used for widening calcified roots in the coronal third.

P4D
PD4D



Used to remove subgingival calculus.

P5



Curved to the left, used to remove subgingival calculus.

P6L
PD6L



Curved to the right, used to remove subgingival calculus.

P6R
PD6R

EMS SAT SIR NSK

ENDODONTICS



120° angle holder for U files with a shank diameter of 0.57mm, this tip is usually used for anterior root canal cleaning.

E1 ES1 EN1
ED1



95° angle holder for U files with a shank diameter of 0.57mm, this tip is usually used for molar root canal cleaning.

E2 ES2 EN2
ED2



With cooling spray disabled, this tip is used for lateral thermal condensation of gutta percha.

E3 EN3
ED3



Diamond coated, used to remove calcification and bad filling material from pulp cavity and overhang pulp chamber.

E3D EN3D
ED3D



Quickly removes obstructions and broken instruments in the mid and upper thirds of large and straight canals.

E4 ES4 EN4
ED4



Diamond coated, for retreatments or extremely hard materials.

E4D EN4D
ED4D

EMS SAT NSK



Useful for quickly eliminating calculus and temporary cements.



Diamond Coated, for enhanced cutting and abrasion. Especially useful for removing very hard materials on the canal wall.

E5

EN5

ESD

EN5D

ED5

ED5D



Bur holder used to expand the root canal and grind teeth.



Bur holder used to expand the root canal and grind teeth.

E8

EN8

E9

EN9

ED8

ED9



Used for gentle treatment of the root canal of anterior teeth.



Diamond coated, used for efficient apical polishing of anterior teeth.

E10

EN10

E10D

EN10D

ED10

ED10D



Used for gentle treatment in the root canal of anterior teeth.



Diamond Coated, with enhanced cutting and abrasive power; especially useful for removing very hard materials on the canal wall and lateral cutting.

E11

EN11

E11D

EN11D

ED11

ED11D

EMS SAT SIR MEC NSK



Quickly removes obstructions and broken instruments in the mid and upper thirds of large and straight canals.



Diamond coated for retreatments or extremely hard materials.

E14

E14D

ED14

ED14D



Useful for quickly eliminating calculus and temporary cements.



Diamond coated, for enhanced cutting and abrasion. Especially useful for removing very hard materials on the canal wall.

E15

E15D


ED15

ED15D

CROWN REMOVAL



Used to remove bridges, crowns and posts by fracturing the cement.

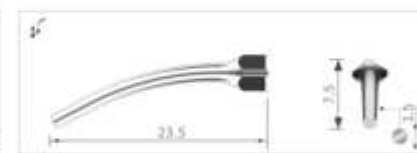

Used to remove bridges, crowns and posts by fracturing the cement.

G7

GS7

GN7

G8

GD7

CAVITY PREPARATION



Ball diamond tip, for preparation of occlusal surfaces and cervical ridges.

SB1

SBD1

SBN1




Semi-spherical diamond tip, for finishing cervical interproximal areas without causing lesions to adjacent teeth.

SB2

SBD2

SBN2

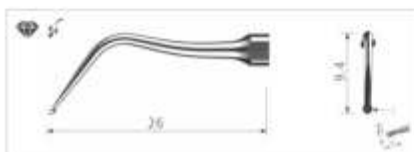


Semi-spherical diamond tip, for finishing cervical interproximal areas without causing lesions to adjacent teeth.

SB3

SBD3

SBN3

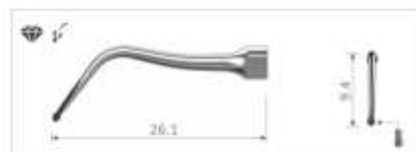

Semi-spherical diamond tip set 45° to the left. Allows the dentist to remove carious areas of a tooth without damaging the adjacent teeth.

SBL

SBDL

SBNL



Semi-spherical diamond tip set 45° to the right. Allows the dentist to remove carious areas of a tooth without damaging the adjacent teeth.

SBR

SBDR

SBNR

IMPLANT CLEANING



Implant Cleaning. PEEK tip for safe and gentle cleaning of implants and restorations. ICP SET contains 1 ICP WREN, 1 ICP and 2 peeks.

eICP SET

sICP SET

SirICP SET

nICP SET

ICPpeek


ICPpeek

ICPpeek

ICPpeek



Replacement. PEEK tip pack of 10

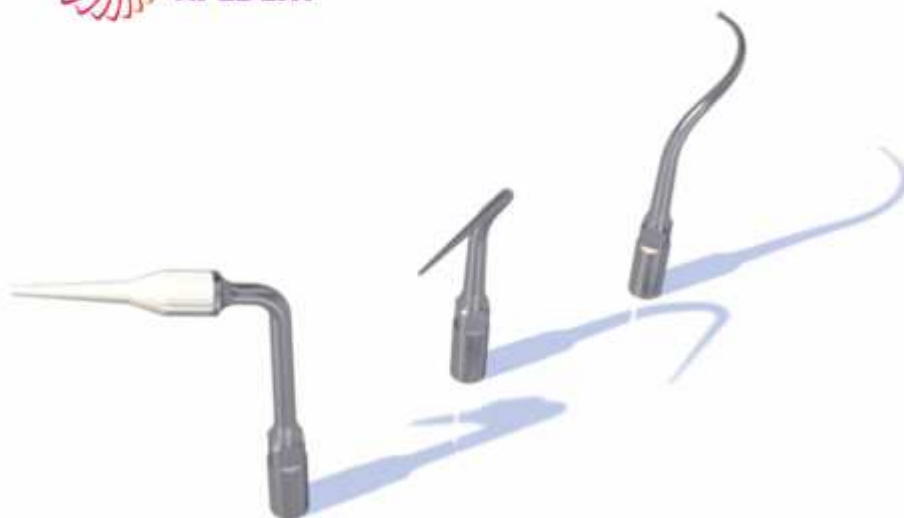

EMS
STYLE
TIPS

SCALING & PERIODONTICS



Universal. Accurate removal of solid supragingival calculus in all quadrants.



Perio. The P instrument is well suited for the precise removal of stubborn subgingival and supragingival calculus and concretions in all quadrants.



PerioSlim. The fast and minimally-invasive PS instrument is ideal for calculus removal in deep periodontal pockets in all quadrants.



Perio Cleaning. For cleaning and disinfecting periodontal pockets, it is especially suitable for follow-up and recall visits for periodontal maintenance patients.



IMPLANT CLEANING

eICP-SET



Implant Cleaning. PEEK tip for safe and gentle cleaning of implants and restorations. eICP SET contains 1 ICP WREN, 1 eICP and 2 PEEKs.

ICPpeek



10X

Replacement. PEEK tip pack of 10

TITANIUM PERIODONTICS

P50L



Curved to the left used for subgingival scaling, periodontal pocket probing and root planing.

P50R



Curved to the right used for subgingival scaling, periodontal pocket probing and root planing.

P52



Used for subgingival scaling, periodontal pocket probing and root planing.

P53L



Curved to the left used for subgingival scaling, periodontal pocket probing and root planing.

P53R



Curved to the right used for subgingival scaling, periodontal pocket probing and root planing.

P54L



Curved to the left used for subgingival scaling, periodontal pocket probing and root planing.

P54R



Curved to the right used for subgingival scaling, periodontal pocket probing and root planing.

P59



Used for subgingival scaling, periodontal pocket probing and root planing.



SAT


SCALING



Universal. Simple supragingival scaling and tartar removal without damage to the enamel.



Calculus. Apply the flat part to the tooth to remove major calculus deposits.

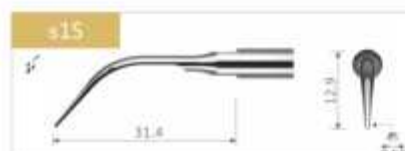

Slim. Supra and sub-gingival scaling. Particularly effective for interproximal spaces.



SATELEC
STYLE
TIPS

IMPLANT CLEANING



Implant Cleaning. PEEK tip for safe and gentle cleaning of implants and restorations. ICP SET contains 1 ICP WREN, 1 ICP and 2 peeks.



Replacement. PEEK tip pack of 10


PERIODONTICS

s10X



Interproximal. Its flat active part makes it suitable for the interproximal spaces and supra-gingival scaling. Its anatomical design enables fast and efficient use.

s10P

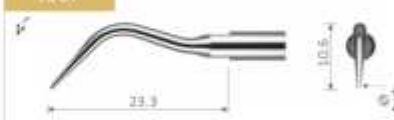

Shallow. Slim round instrument suitable for scaling shallow pockets less than 2-3mm deep.

s10Z



Sub-Gingival. Medium pockets (< 4mm). Removal of biofilm and soft deposits, enables observation of pocket depth using the 3mm scale.

sH3



Anterior periodontics. Enables treatment of the incisor-canine block. The guide edge is oriented parallel to the pocket. Enables work in the periodontal pocket without injury to the ligament.


ENDODONTIC

sET40D



Diamond. A diamond-coated steel tip for retreatment of very hard material in the middle third. Length 40mm, 30 µm diamonds, taper 4%.

sH4L

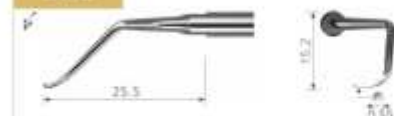


Premolar and molar, left-curve. First instrument for treating all the surfaces and the furcations.
Maxillary: buccal and distal surfaces of sector 2, pivots at 13, then the buccal and mesial surfaces of sector 1.
Mandibular: buccal and distal surfaces of sector 4, pivots at 43, then lingual and mesial surfaces of sector 3.

sH4R



Premolar and molar, right-curve. second instrument, use after H4L.
Maxillary: palatine and mesial surfaces of sector 2, pivots at 13, then buccal and distal surfaces of sector 1.
Mandibular: lingual and mesial surfaces of sector 4, pivots at 43, then buccal and distal surfaces of sector 3.

sTK1L

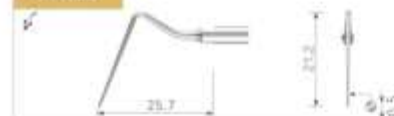

Extended Probe. Ideal for examining and maintenance of medium to deep pockets (> 4mm). A diagnosis aid used during the debridement and irrigation of pockets.

sTK1S




Probe. Recommended for examining shallow and medium pockets (< 4mm) and for the maintenance of simple cases. sTK1S tip is graduated every 3mm.

sTK2L



Left-Oriented Maintenance. Recommended for the maintenance of moderate to deep pockets and furcations. Equivalent to the Nabers probe.

sTK2R

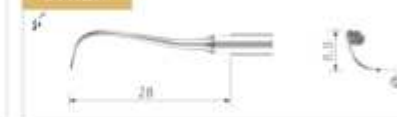

Right-Oriented Maintenance. Partner to the sTK2L. Recommended for the maintenance of moderate to deep pockets and furcations. Equivalent to the Nabers probe.



NSK

AIR SCALING



NSK
STYLE
TIPS

nS1



Universal. Removal of marginal-gingival and supra-gingival calculus.

nS2



Perio. Removal of supra-gingival calculus.

nS3



Sickle. Removal of supra-gingival calculus.

NK4



Universal slim tip.

NK5



A root cleaning tip, straight application.

NK6



A root cleaning tip, with a left curve.

NK7



A root cleaning tip, with a right curve.

nSICP SET



Implant Cleaning. PEEK tip for safe and gentle cleaning of implants and restorations.
nSICP SET contains 1 ICP WREN, 1 nSICP and 2 PEEKs.



SIRONA STYLE TIPS



SIRONA

SCALING

sr1L



General. A flat tipped instrument for initial lingual and buccal scaling.

sr2L



Versatile. A versatile instrument used for lingual, buccal and approximal scaling.

sr3L



Interdental. A universal instrument which can be used for subgingival scaling at interdental sites.

sr4L



Universal. A fine instrument, ideal for supragingival and subgingival tooth surfaces and interdental areas.

IMPLANT CLEANING

sirICP SET



Implant Cleaning. PEEK tip for safe and gentle cleaning of implants and restorations.
sirICP SET contains 1 ICP WREN, 1 sirICP and 2 peeks.

ICPpeek



Replacement. PEEK tip pack of 10



KAVO STYLE TIPS

KAVO

SONICflex

k#5



General. For calculus removal from all tooth surfaces.

k#6



Fine. Fine instrument for smaller areas and those which are difficult to access.

k#7



Subgingival. Instrument used for subgingival scaling.

k#8



Subgingival. For the removal of subgingival calculus.

k#60



A long straight tip with rounded surfaces, the k#60 is used for gentle pocket cleansing during periodontal treatment. It is effective cleansing in pockets up to 9 mm deep. It is used in the anterior region, smooth surfaces.

GK2



Used for supragingival, subgingival and interdental calculus removal.

GK5



Used for supragingival, subgingival and interdental calculus removal.

GK6



Used for supragingival, subgingival and interdental calculus removal.

GK7



Used for supragingival, subgingival and interdental calculus removal.

EK8



Bur holder, used to expand the root canal.

kICP SET



Implant Cleaning. PEEK tip for safe and gentle cleaning of implants and restorations.
kICP SET contains 1 ICP WREN, 1 kICP and 2 peeks.

ICPpeek



Replacement. PEEK tip pack of 10

PIEZOlux & SONOsoft

k#9 (GC2)



Used to remove calculus and bacterial plaque from supragingival interdental calculus.

k#10 (GC1)



Used to remove calculus and bacterial plaque from supragingival interdental calculus.

k#11 (PC1)



Used to remove subgingival calculus.

SONICflex Quick 2008

k5A



General. For calculus removal from all tooth surfaces.

k6A



Fine. Fine instrument for smaller areas and those which are difficult to access.

k7A



Subgingival. Instrument used for subgingival scaling.

k8A



Subgingival. For the removal of subgingival calculus.

k60A



A long straight tip with rounded surfaces, the k#60 is used for gentle pocket cleansing during periodontal treatment. It is effective cleansing in pockets up to 9 mm deep. It is used in the anterior region, smooth surfaces.

KICPA SET



Implant Cleaning. PEEK tip for safe and gentle cleaning of implants and restorations. KICPA SET contains 1 ICP WREN, 1 KICPA and 2 peeks.

ICP peek



Replacement. PEEK tip pack of 10

PIEZOsoft & PiezoLED

k201



Universal. Accurate removal of solid supragingival calculus in all quadrants.

k202



Perio. The P instrument is well suited for the precise removal of stubborn subgingival and supragingival calculus and concretions in all quadrants.

k203



PerioSlim. The fast and minimally-invasive PS instrument is ideal for calculus removal in deep periodontal pockets in all quadrants.



AMDENT STYLE TIPS

AMDENT

SCALING



Universal tip for subgingival scaling, furcations, supragingival fine scaling and spot removal.



Primerily for lingual and buccal subgingival scaling and furcations. Also used for supragingival fine scaling and spot removal.



Universal tip, can be used for lingual and buccal supragingival scaling.